


# PERFORMALIFT SERIES

Our most comfortable table just got better! The Performalift has always offered the ultimate in comfort with the ability to tilt the table up at the head end to relieve sinus pressure or tilt it in the opposite direction with a salon or backrest top to create a zero gravity chair position. Now standard with a removable and replaceable padded top!



**New!**  
One Piece Removable Pad

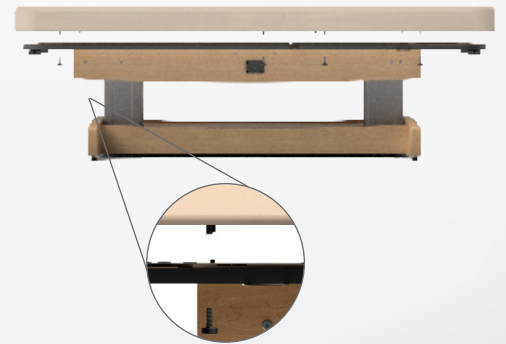

Padding and fabric on tables never last as long as the bases, and occasionally they get damaged; in the past you had to replace the whole table. Our new design features finger knobs that can be removed in minutes, allowing you to remove the foam and fabric and install a new one without any special tools so your table will be serviceable for much longer.



## FEATURES




### Zero Gravity

The unique twin tower lift system creates a Zero Gravity chair capability insuring the ultimate in comfort. The tilting function is helpful for many treatments and offers greater comfort for people who experience sinus pressure when lying flat.



### Programmable Hand Control

Now standard, this control has two user programmable memory positions and the ability to change two settings at once with a single button. The simple interface makes positioning clients fast and the new electronics are quiet and smooth in operation.

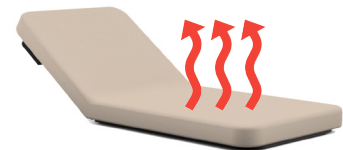

### Integrated Heated Top (Optional)

Our new system is imbedded in the foam, creating a more even heat, and has three settings. Designed originally to be used in car seats and hospital beds, the heating element automatically shuts off after 10 hours providing peace of mind for safety and energy conservation.

made in the US



**OAKWORKS®** Spa

# SPECIFICATIONS

- New! Removable and replaceable padding
- Tilttable top for greater comfort
- Programmable memory positions
- AeroCel™ padding for ultimate comfort
- Large height range - 20" - 34"

## TABLE SPECIFICATIONS

Lift System	Electric, Twin Tower
Controls	Foot Control + Hand Control with memory positions
Lifting Capacity	550 lbs. (249.5 kg)
Fabric	TerraTouch - Medical Grade, ISO 10993
Finish	Water Based Stain or Lacquer
Padding	4" (10.2 cm) or 6" (15.2 cm) AeroCel™ (Optional)
Height Range 4" (10.2 cm) pad	20" - 34" (51 cm - 86 cm)
Tilt	+/- 18°
Width	29" (73.6 cm), 31" (78.7 cm), 33" (83.8 cm)
Length	73" (185.4 cm)
Options (must be ordered with table)	Integrated Heated Top, ABC System, Casters
Warranty	3 years parts, 2 years labor
Safety Listings	cETLus listed 3034177
Accessories	Adjustable Side Arm Rests, Boiance Shelf Support, Boiance Face Rest Cushion, Quick Lock Face Rest

## STAINS



## LACQUERS



## TERRATOUGH COLORS

Meets the requirements of ISO 10993-1:2009 Biological Evaluation of Medical Tables & passes testing for Cytotoxicity, Skin Irritation & Sensitization as well as CA TB1 17 Fire Retardance testing. Please refer to the product cleaning instructions for a list of approved cleaners and disinfectants.




## ABC SYSTEM (Optional)

Oakworks Patented ABC system allows you to create the perfect level of support for each client, taking pressure off their breasts and tension out of their back. The system also raises the pad above level, making it a great head support for face up positions, and can be used for stretches also.



Incredible breast comfort  
Lower to lessen pressure

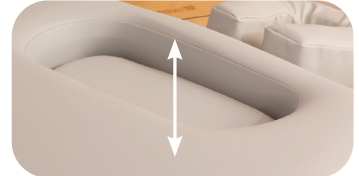

Electric adjustments  
Large adjustment range



Head support  
Built in head support



Pectoral muscle stretch  
Advanced positioning

## ADJUSTABLE SIDE ARMREST

These versatile and comfortable removable arm rests feature a wide range of positions. They always stay level, regardless of the backrest angle. They offer comfortable positioning and can be moved completely out of the way to allow guests to get on and off the table easily.



## BOIANCE SHELF SUPPORT

The shelf and bracket used together allow you to create the perfect place for clients arms when working on their back. The shelf can also be used alone as a head support for facials and other aesthetic treatments. The shelf is designed to allow therapists to relax their shoulders while working on clients, minimizing stress and injuries.



717.235-6807 [www.spatables.com](http://www.spatables.com) email: [customerservice@oakworks.com](mailto:customerservice@oakworks.com)